

MyPipetman®

SELECT

PRECISION THAT ADAPTS TO YOUR LAB

Ergonomic, safe, and reliable, MyPIPETMAN® Select offers volume setting choice and optional text personalization on the handle for an optimal use in various protocols.

Set and Secure the Volume your Preferred Way

The unique Trilock™ volume-locking system enables you to choose different positions when setting the volume:

- Free volume setting position for quick and effortless volume adjustment
- Light ratchet volume setting for precise volume adjustment
- Volume locked position for secure pipetting

With the Trilock system, the volume is adjusted simply by turning the push button of the pipette.



Ergonomic Design

Lightweight and well-balanced, MyPIPETMAN Select has an elegant design with a textured effect for a better grip.

With no compromise on quality and performance, the built-in ejection system provides reduced tip ejection forces. This patented ejection system reinforces the pipette's ergonomics for convenient use throughout the day.

Personalize your Pipette for Optimal Use

MyPIPETMAN Select can have their handles printed on demand. Identify your pipette to avoid any risk of contamination by dedicating it to a specific liquid or workflow step and improve your lab Lean Management process.



MyPIPETMAN Select can be personalized with a special text. Part numbers are finishing with SP (Select Personalized) instead of S (Select). Please contact your [Gilson Account Manager](#) for special text printing and SP part numbers.

Extra Customization for Your MyPIPETMAN®.

For identification or more convenient use, this colored window will provide extra customization to your MyPIPETMAN®. Six colors are available, and you can easily replace the colored window at your convenience.



Colored Windows

Model	Part Number
Window 3 digit Dark Blue MyPIPETMAN	FP070101
Window 3 digit White MyPIPETMAN	FP070102
Window 3 digit Blue MyPIPETMAN	FP070103
Window 3 digit Red MyPIPETMAN	FP070104
Window 3 digit Yellow MyPIPETMAN	FP070105
Window 3 digit Green MyPIPETMAN	FP070106

SHOP NOW




WATCH VIDEO

Store Your Pipettes On Our CARROUSEL™ Pipette Stand

The CARROUSEL™ Pipette Stand holds up to seven Gilson single or multichannel pipettes, protecting them from contamination.

Full range of accessories can be found on [page 74](#).

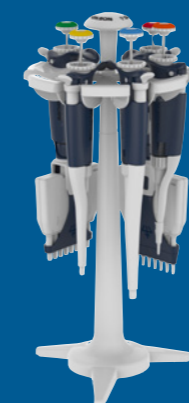

SHOP NOW



MYPIPETMAN® SELECT SINGLE CHANNEL MODELS

Model	PIPETMAN® DIAMOND Tips	Part Number	Volume Range	
P2	D10 DL10	DF10ST DFL10ST	FP10001S	0.2-2 µL
P10	D10 DL10	DF10ST DFL10ST	FP10002S	1-10 µL
P20	D200	DF30ST	FP10003S	2 -20 µL
P100	D200	DF100ST	FP10004S	10-100 µL
P200	D200 D300	DF200ST DF300ST	FP10005S	20-200 µL
P1000	D1000 D1200	DF1000ST DF1200ST	FP10006S	100-1000 µL
P5000	D5000		FP10007S	500-5000 µL
P10mL	D10mL		FP10008S	1000-10000 µL

MYPIPETMAN® SELECT MULTICHANNEL MODELS

P8x10	D10 DL10	DF10ST DFL10ST	FP10009S	0.5-10 µL
P12x10			FP10010S	
P8x20	DL10 D200	DFL10ST DF30ST	FP10011S	2-20 µL
P12x20			FP10012S	
P8x200	D200 D300	DF200ST DF300ST	FP10013S	20-200 µL
P12x200			FP10014S	
P8x300	D200 D300	DF200ST DF300ST	FP10015S	20-300 µL
P12x300			FP10016S	

SHOP NOW

MYPIPETMAN® SELECT KIT

Model	Part Number	Models	PIPETMAN® DIAMOND Tips	Accessories
MyPIPETMAN® Select 4-Pipette Kit	F167390	P2 P20 P200 P1000	DL10 D200 D1000	4 SINGLE pipette holders 4 Plastic ejectors 6 colored windows (1 of each color)

SHOP NOW



Tools and Resources

Specsheet: Pipette Specifications

DOWNLOAD

Guide to Pipetting: Tips from Our Experts

DOWNLOAD

Video: MyPIPETMAN Two-Minute Inspection Video

WATCH